




An ISO 9001:2015, EN 13485:2012, OHSAS 18001 :2007 &  Certified Company

Remote Operated Obstetric OT Table



Remote Operated Obstetric OT Table CELEBRITA

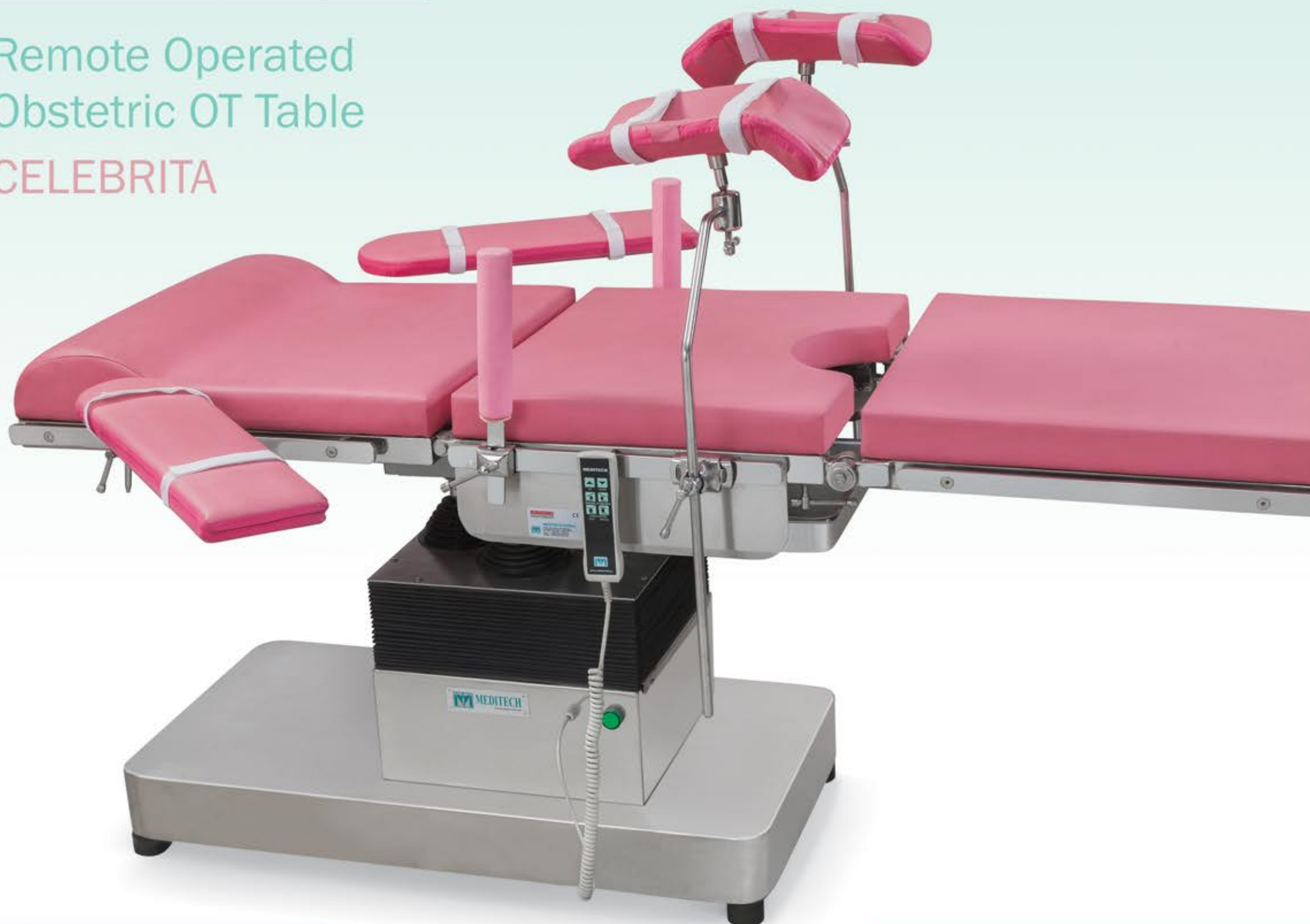

Height adjustment



Trendelenburg / Rev. Trend.



Back section adjustment



Lithotomy Position

Features

- Electromatic drives with noiseless, jerk less and smooth positioning.
- Three section upholstered table top ensures maximum comfort to the patient.
- Table positions i.e. Height adjustment, Trendelenburg, Reverse Trendelenburg, Back Section/Chair Position are controlled by corded handset with graphical icons.
- Detachable leg sections offer perineal access & allow to attach drainage tray.
- Integrated weight compensation by manual adjustment of leg section.
- Large size stainless steel drainage tray.
- Non-corrosive stainless steel 304G construction including table top frames, accessories rail, base cover and column telescopic covers.
- Easy for cleaning & hygiene.
- The powered movements are activated by imported linear actuators.
- Four easy glide, antistatic castors, two of which are swivelling, ensure mobility.
- Pedal operated central braking system secures the table in its place while in use.

Standard accessories

- Lithotomic crutches – 1 pair
- Drainage tray – 1 pc
- Arm board – 1 pair
- Hand grip – 1 pair

TECHNICAL CONSTRAINTS

Top dimension	L 1750 x W 600 mm
Height adjustment with pad	700mm - 950 mm
Trendelenburg / Rev. Trend.	15° / 15°
Back rest (upward / downward)	60° / 0°
Leg rest (upward / downward)	15° / 90°
Power supply	24V DC
Patient weight capacity	150 Kgs (330 lbs)

Major Clients



Presence across major parts of the World



MEDITECH INDIA®

487/67, Glass Street, National Market, Peeragarhi, New Delhi-87, India.

Tel. : +91-11-40392484

Mob. : +91- 9812266676, 9812035477

Email : info@meditech-india.com marketinguk@meditech-india.com

Web : www.meditech-india.com www.surgical-lights.com

Skype : exportmeditech

Authorised Distributor

Design & specifications are subject to change without any further notice.

Publishing of any material in full or half of this catalogue is a punishable offence, as the products shown in the catalogue are the sole property of MEDITECH INDIA.